




Abbott

ENABLE TIMELY DIFFERENTIAL DIAGNOSIS AND TREATMENT OF COVID-19 AND FLU A&B USING A SINGLE SWAB


PANBIO™ COVID-19/FLU A&B RAPID PANEL (NASOPHARYNGEAL)

- For patients presenting symptoms or suspected of current COVID-19, Flu A and/or B infection
- Only 1 nasopharyngeal swab for both COVID-19 and Flu A&B
- Simplified workflow with pre-filled buffer tube
- No instruments needed to interpret results
- Results in 15 minutes

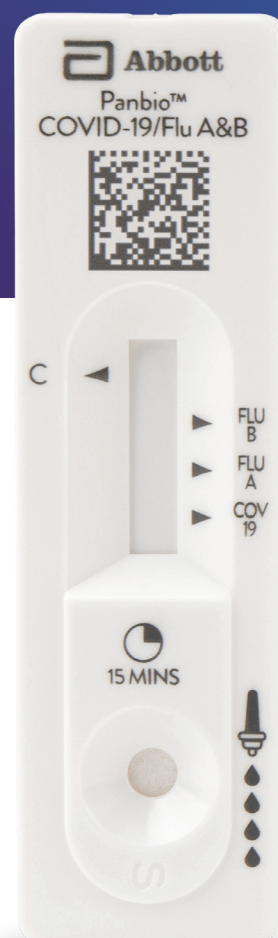

Abbott is the trusted brand used in more **100 countries**.


DESIGNED FOR EASY TESTING

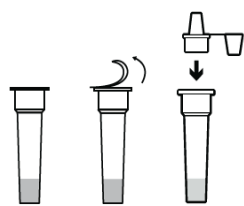

- PRE-FILLED BUFFER TUBE WITH INTEGRATED DROPPER CAP. NO MANUAL FILLING NEEDED




- ONLY ONE SPECIMEN COLLECTING SWAB NEEDED

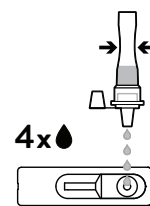

- SIMPLE TEST CASSETTE DESIGN WITH ONE SAMPLE WELL

TEST PROCEDURE

Consult Instructions for Use for complete procedure.

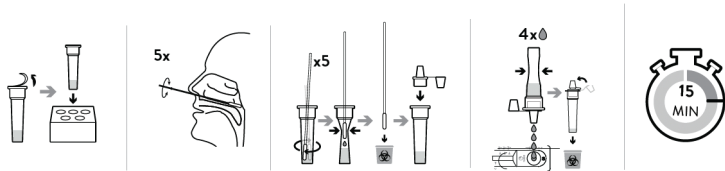

SPECIFICATIONS

- **TEST TIME:** 15 MINUTES
- **STORAGE:** 2°C–30°C
- **CE MARK**
- **SAMPLE TYPE:** Nasopharyngeal (NP) swab

PERFORMANCE

	COVID-19	COVID-19 (Ct≤30)	FLU A	FLU B
SENSITIVITY	88.0%	94.4%	92.0%	100%
SPECIFICITY	100%	100%	100%	100%
95% CI (Sens.)	80% - 93.6%	87.5% - 98.2%	80.8% - 97.8%	92.9% - 100%
95% CI (Spec.)	98.2% - 100%	98.2% - 100%	98.2% - 100%	98.2% - 100%

INTENDED USE: The Panbio™ COVID-19/Flu A&B Rapid Panel (NP) is an *in vitro* diagnostic rapid test for the qualitative detection of SARS-CoV-2 and Influenza A&B antigen (Ag) in human nasopharyngeal swab specimens from individuals, who meet COVID-19 and Influenza clinical and/or epidemiological criteria. Panbio™ COVID-19/Flu A&B Rapid Panel (NP) is for professional use and is intended to be used as an aid in the diagnosis of SARS-CoV-2 and Influenza infection. The product may be used in laboratory and non-laboratory environment.

The test provides preliminary test results. Negative results don't preclude SARS-CoV-2 and Influenza A&B infection and they cannot be used as the sole basis for treatment or other management decisions. Negative results should be combined with clinical observations, patient history, and epidemiological information.

ORDER INFORMATION

- **Panbio™ COVID-19/Flu A&B Rapid Panel (Nasopharyngeal)**
- **CATALOG NUMBER:** 62FK10; 10 tests/kit (control kit sold separately; catalog number: 62FK11C)
- **CONTENTS:** 10 Test Devices, 10 Extraction Buffer Tubes, 10 Nozzle Caps, 10 Sterilized Nasopharyngeal Swabs For Specimen Collection, 1 Tube Rack, 1 Instructions For Use (IFU) with Quick Reference Guide (QRG)

CONTACT YOUR LOCAL REPRESENTATIVE TODAY.

globalpointofcare.abbott

Product not available in all countries.

